



ENTERAL  
Enteral feeding pumps and accessories



easyfit & easyMoov6

## Enteral feeding pumps

Easy to use and easy to set up, **easy**lit**** & **easy**moov6**** enteral feeding pumps are intuitive and adapted to stationary and ambulatory needs in hospital and homecare.

### easy**lit**

A simple pump, designed for stationary use



#### Intuitive

Short learning curve

#### Powerful

A continuous rotary peristaltic rotor helps delivering high viscosity food

#### Quiet

A continuous rotary peristaltic rotor allows the pumps to be quiet. This is especially appreciated for night feeding

#### Easy to clean

Pumps rounded shapes allow easy decontamination

- Weight: 1400 g
- Autonomy: 7 hours at 125 ml/h
- Accuracy: 10%
- Flow rate range: 1 to 600 ml/h
- Volume range: 1 to 9,999 ml
- Automatic prime: 1 min 35 sec
- Memory: 24 hours
- Visual and audible alarms
- Integrated pole clamp

### easy**moov6**

An intuitive pump, designed for ambulatory use




- Weight: 650 g
- Autonomy: 20 hours à 125 ml/h
- Accuracy: 5%
- Flow rate range: 1 to 400 ml/h
- Volume range: 1 to 9,999 ml
- Automatic prime: 1 min 45 sec
- Memory: 24 hours
- Visual and audible alarms
- Multiple languages available
- **easysoft**

Reference	Designation	Qty/box
0VEPL4Xxx <sup>(1)</sup>	easy <b>lit</b> 4 stationary pump: programmable flow rate	1
0VEPL6Xxx <sup>(1)</sup>	easy <b>lit</b> 6 stationary pump: programmable flow rate and volume	
0VEPM6Xxx <sup>(1)</sup>	easy <b>moov6</b> ambulatory pump: programmable flow rate and volume, nutrition mode	

<sup>(1)</sup> Xxx : required plug and language version to specify

## Giving sets

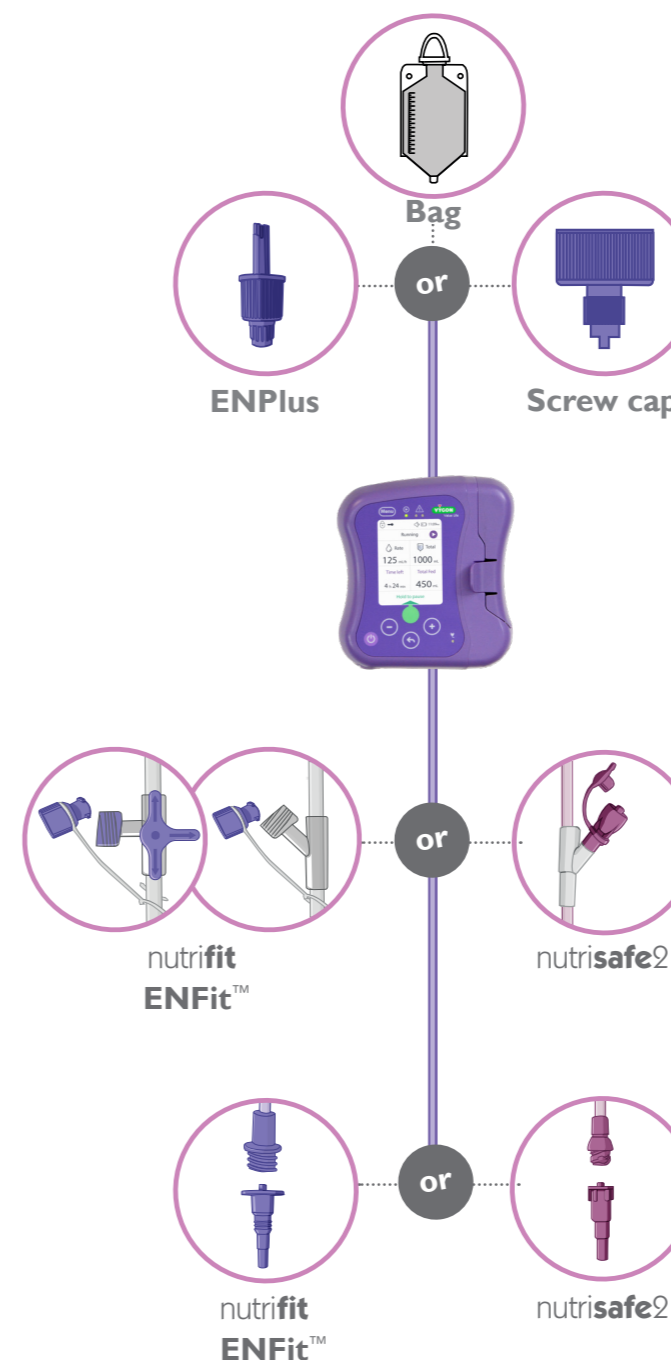

### Proximal connection

ENPlus connection mitigates risks of misconnection with parenteral solution bags.

Enteral nutrition bags are available in 500 ml, 1200 ml and 1500 ml.

The screw cap connection allows a direct connection to bottles.


### Lateral connection

Lateral connection of giving sets is available with the safe connection ENFit™ or with Vygon's safe system **nutrisafe2**.<sup>(2)</sup>


### Distal connection

Distal connection of giving sets is available with the safe connection ENFit™ or with Vygon's safe system **nutrisafe2**.<sup>(2)</sup>

Giving sets are delivered with a christmas tree connector screwed on the distal connection.

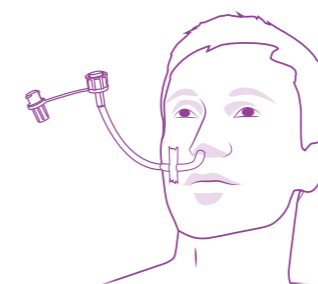

<sup>(2)</sup> Vygon recommends **nutrifit** (ENFit™) for paediatric and adult use and **nutrisafe2** for neonatal use.

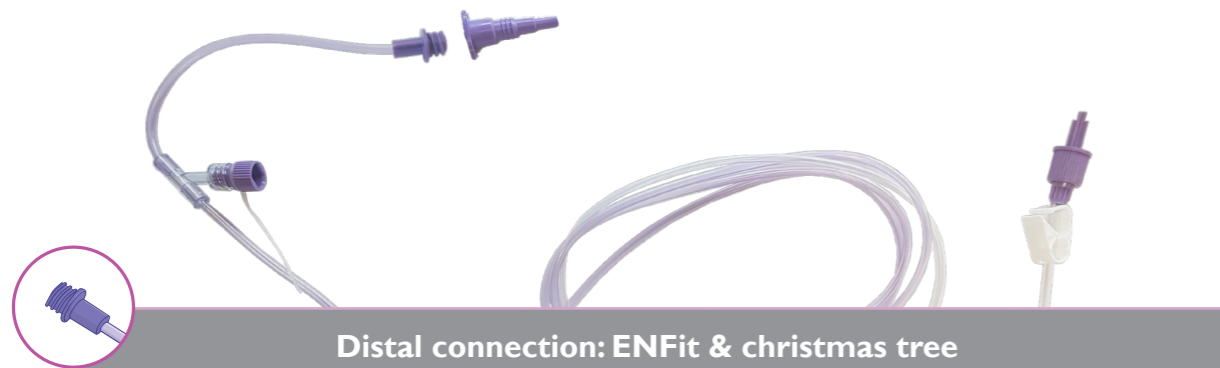

Distal connection: ENFit & christmas tree

Pump Proximal connection	easylit <sup>(3)</sup>		easymoov6	
	Reference	Lateral connection	Reference	Lateral connection
	0VENEJ34S	ENFit	0VENEP14S	ENFit
			0VENEP08S with valve	ENFit
			0VENEP14R	3-way stopcock
			0VENEP11S	-
	0VENEJ34SCS	ENFit	0VENEP14SCS	ENFit
			0VENEP14RCS	3-way stopcock
			0VENEP11SCS	-
	0VENP2GBJ34S 1200 ml	ENFit	0VENP1P14S 500 ml	ENFit
			0VENP2P14S 1200 ml	ENFit
	0VENEJ24S	ENFit	0VENEJ24S	ENFit
	0VENEJ24SH	ENFit	0VENEJ24SH	ENFit
	0VENEPJ24SHCS 1200 ml	ENFit	0VENEP2P14SHCS 1200 ml	ENFit
	-	ENFit	0VENP2P14SH 1200 ml x 2	ENFit



Distal connection: nutrisafe2 & christmas tree

Pump Proximal connection	easylit <sup>(3)</sup>		easymoov6	
	Reference	Lateral connection	Reference	Lateral connection
	0VENEJ34V	nutrisafe2	0VENEP14V	nutrisafe2
			0VENEP08V with valve	
	0VENEJ34VCS	nutrisafe2	0VENEP14VCS	nutrisafe2
	0VENP1GBJ34V 500 ml	nutrisafe2	0VENP1P14S 500 ml	nutrisafe2
	0VENP2GBJ34V 1200 ml			
	0VENEJ24V	nutrisafe2	-	-

<sup>(3)</sup> easylit giving sets are also compatible with the following pumps: Jumbo, Nestlé P7000, Compat Standard, Flocare 800.

## Bags

All bags are fitted with:

- A connection compatible with ENPlus
- An ice compartment for cooling the enteral formula
- A transparent front side, which allows to write the patient and formula's information.

Reference	Designation	Qty/box
0VENP1GT	500 ml bag	50
0VENP2GT	1200 ml bag	
0VENP3GT	1500 ml bag	



## Giving set accessories

A broad range of accessories is available to customize enteral giving sets according to clinical needs.

Reference	Designation	Qty/box
0VEACC015	Universal bag adaptor	300
0VEACC023	Clamp	50
0VEACC126	Upper double line adaptor with hydration line and ENPlus spike adaptor	100
0VEACC127	Upper double line adaptor with ENPlus spike adaptor	100
0VEACC406	ENFit™ / christmas tree adaptor	50
0VEACC416	ENFit™ / christmas tree adaptor with cap	50
0VEACC500	ENFit™ / funnel adaptor	50
0VEACC510	ENFit™ / nutrisafe2 adaptor	50
368.72	nutrisafe2 / ENFit adaptor	100
368.12	Taper connector with nutrisafe2 male connection	30
368.02	nutrisafe2 / feed set adaptor	30
829.01	Cap for nutrisafe2 feeding tubes	100

## Pumps accessories



Reference	Designation	Qty/box
0VEPM6BAG	easy <b>moov6</b> backpack	1
0VESUPPOMP	easy <b>moov6</b> pump table stand	1
1MWPM6CHA <sup>(4)</sup>	easy <b>moov6</b> AC adaptor	1
1MWPM6ACP	easy <b>moov6</b> pole clamp	1
1MWPLCHA02 <sup>(4)</sup>	easy <b>lit</b> AC adaptor	1

<sup>(4)</sup> : required plug version to specify

## easysoft, exclusively for easymoov6

easy**soft** is a software to help in the management of the easy**moov6** pump allowing:

- Access to the 500-events history calendar of the pump
- Customization of some pump parameters
- Help in the preventive maintenance of the pump

easy <b>moov6</b> pump			
Reference	Designation	Additional information	Qty/box
1MWEASYSOFT	easy <b>soft</b>	Delivered as USB key with USB-Micro USB cable	1



**Important:** easy**soft** is available in French and English **only**.

 OBSTETRICS NEONATOLOGY ENTERAL

**For further information, please contact: [questions@vygon.com](mailto:questions@vygon.com)**

The specifications shown in this leaflet are for information only  
and are not, under any circumstances, of a contractual nature.

---

Vygon – 5, rue Adeline • 95440 ECOUEN • FRANCE  
Reception: +33 (0)1.39.92.63.63 – Service clients France: +33 (0)1.39.92.63.81  
Export customer service: +33 (0)1.39.92.64.15  
Fax.: +33 (0)1.39.92.64.44 • [www.vygon.com](http://www.vygon.com)



[www.vygon.com](http://www.vygon.com)